Appl. No.: Unknown

Atty Docket: ABIOS.042A

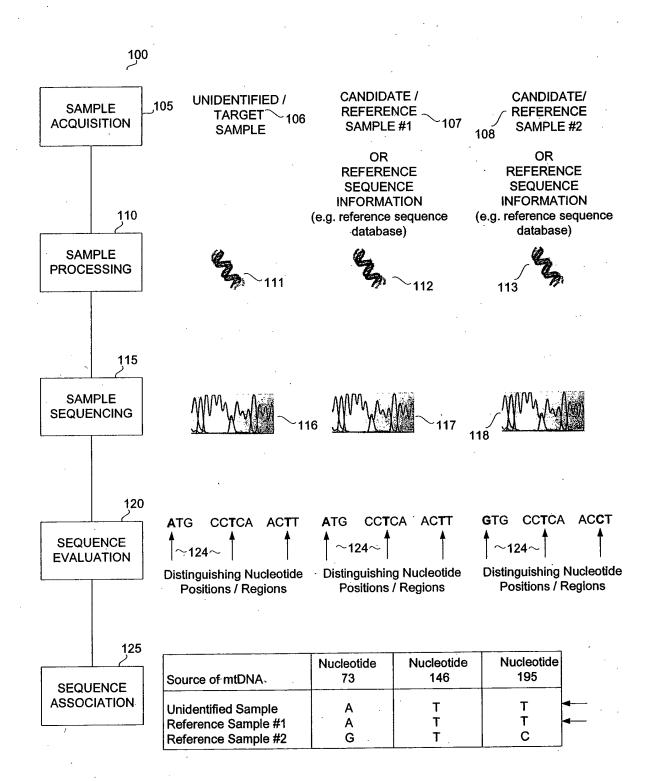


FIGURE 1

Appl. No.: Unknown

Atty Docket: ABIOS.042A


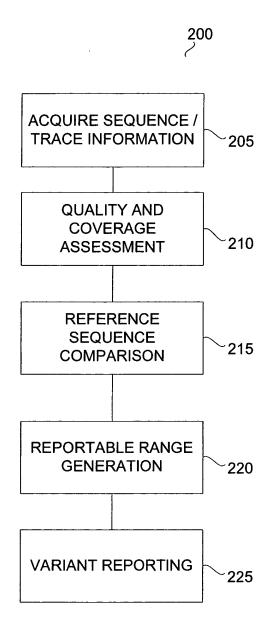


FIGURE 2

Appl. No.: Unknown

Atty Docket: ABIOS.042A

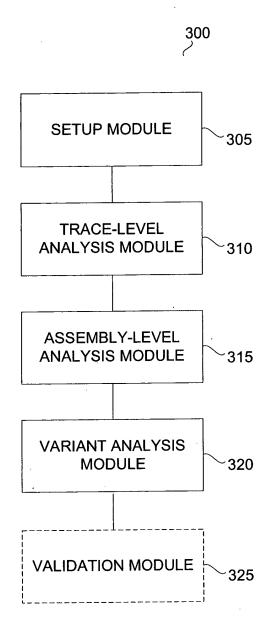

235

FIGURE 3

Atty Docket: ABIOS.042A

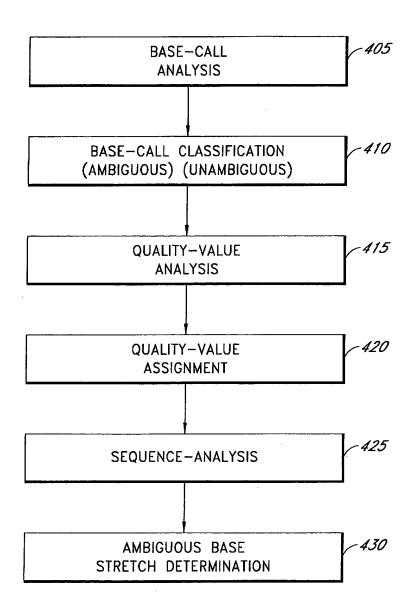


FIG. 4A

Appl. No.: Unknown Atty Docket: ABIOS.042A

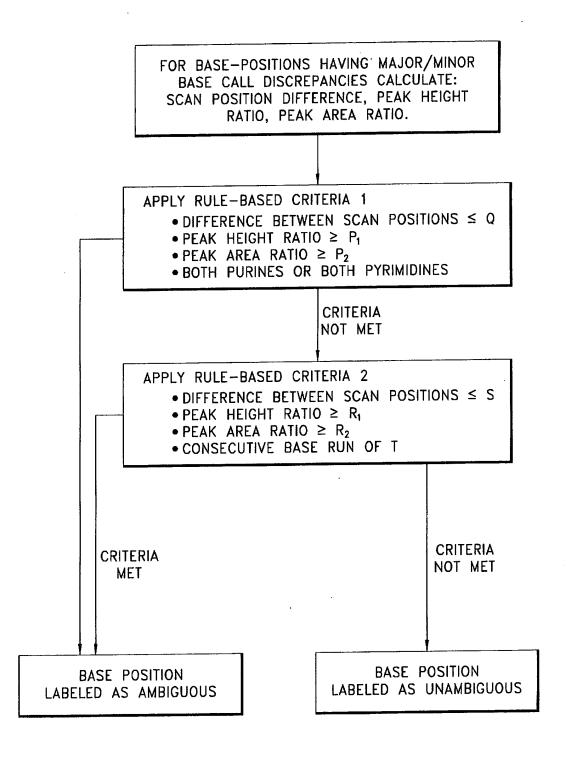


FIG. 4B

Stockwell et al.

Appl. No.: Unknown Atty Docket: ABIOS.042A

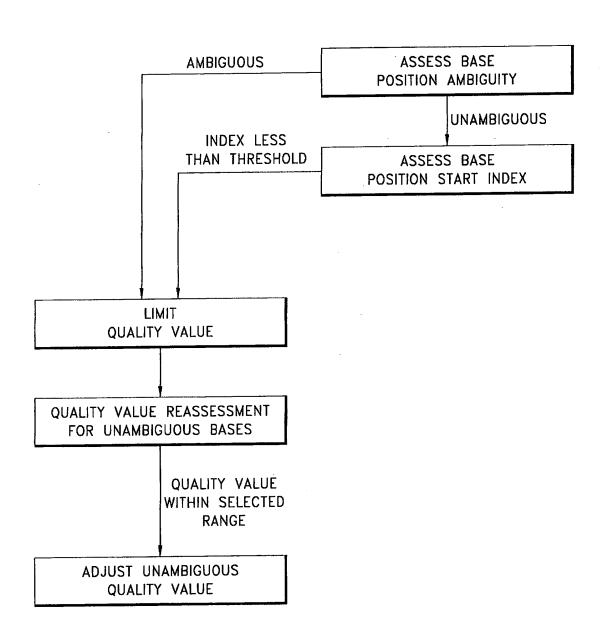


FIG. 4C

Appl. No.: Unknown Atty Docket: ABIOS.042A

Table 1 System components and tunable parameters

Drogram	Function performed		
Program BlastParse.pl	Parse BLAST output		
mark_substitution_heteroplasmy.pl	Mark ambiguous base calls		
extract SE consensus.pl	Extract a consensus sequence from a		
extract_ot_conscrists.pr	phrap polyphred run, replacing low quality		
	bases with '?', and polyphred rank 1		
	bases with 'N'.		
compute coverage.pl	Calculates overall, forward strand, and		
	reverse strand coverage for all bases in		
	the consensus sequence		
count hv1 deletes.pl	Count deletions in HV1 relative to the		
	rCRS		
border_index.pl	Compute start and end positions of HV1		
	and HV2 regions.		
fix_mitotype_reporting_range.pl	Output a mtDNA profile based on input		
	variants list and reportable range. Reject reads that do not align		
find_bad_traces_from_blast_report.pl	110,000		
	appropriately to the rCRS Adds template name, template type, and		
determineReadTypes.pl	primer type to phred output files.		
Odelle	Align mtDNA profile to rCRS and report		
seq2delta_vs	variants according to nomenclature.		
calculate_coverage_mitotype.pl	Mask mtDNA consensus sequence based		
Calculate_coverage_mitotypo.pr	on required coverage and calculated		
	coverage		
count hv1 inserts.pl	Count insertions in HV1 relative to the		
County (_interiorp)	rCRS		
generate_hv_mask_fasta_files.pl	Extract the HV1 portion and HV2 portion		
	of the mtDNA consensus sequence,		
	based on the computed HV1 and HV2		
	regions.		
flip_fasta.pl	Invert an X masked sequence output by		
	cross_match		

MITOCHONDRIAL DNA AUTOSCORING SYSTEM Stockwell et al.

Appl. No.: Unknown

Atty Docket: ABIOS.042A

Table 2

	Assessment 1	Assessment 2
Total True Positives: N _{automated} = N _{manual}	16	10
Total True Negatives: A _{automated} = A _{manual}	77,358	110,354
Total False Positives: N _{automated} = A _{manual}	95	807
Total False Negatives: A _{automated} = N _{manual}	6 ^a	14
Total Incorrect	0	0
Sensitivity: TP/(TP+FN)	72.73%	41.67%
Specificity: TN/(TN+FP)	99.88%	99.27%
Positive Predictive Value: TP/(TP+FP)	14.41%	1.22%
Negative Predictive Value: TN/(TN+FN)	99.99%	99.99%
,		

^a Predominate base correctly called

Differences take into account the consensus of the two analysts